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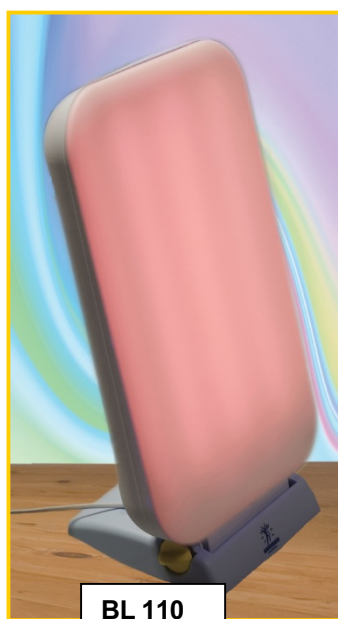
Dear customers,

We are pleased to welcome you as a new customer to our product. You have made the right decision. You have purchased a product which will improve your quality of life for many years.

 Before you use the instrument, please read carefully this user manual.

Effect of light therapy

With the beginning of the age of 25 the production of collagen decreases. This leads to an accelerated skin aging, wrinkles, cellulite, loosing skin and uneven skin tone. Collagen is one of the major proteins of the body and collagen fibers are the most common components of connective tissue.

The **collagen light therapy units** use special wavelengths which achieve exactly the skin layers in which the connective tissue fibers are produced. The collagen production can be stimulated and the collagen in the connective tissues can be increased. So the application of the light therapy can lead to an anti-aging effect with a visible rejuvenation of the skin. Wrinkles and fine lines can be refilled. Pores can be refined, the metabolism can be activated and the whole skin can be improved. The connective tissue can be tightened considerably and the elasticity of the skin can be increased significantly.

How to use the DAVITA® light therapy units

The light therapy units can be used for the Anti-Aging-application to improve the collagen level, to reduce wrinkles, fine lines, blemishes and cellulite. The skin elasticity can be increased.

Some basic information before using the light therapy units

Before use, kindly check that the instrument is in functional order and that there is no sign of damage. The user must check that the equipment functions safely and see that it is in proper working condition before being used. Do not position the equipment so that it is difficult to operate the disconnection device.  Phototherapy instruments must only be operated in dry and closed rooms that are free of any explosion risks. The device is not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide. The instruments must be standing free and should not be exposed to sun.

How to carry out light therapy

Place the unit on a flat surface such as a table and connect the plug with a wall socket.

The Collagen-Light therapy unit BeautyLite BL 70 has on the back side a fixture for folding, so that the device has a save stand. Sit down in front of the instrument and keep a distance of approx. 30 cm. The application duration is in this distance for BL 70 approx. 30 minutes and for BL 110 approx. 20 minutes.

Move the tilting switch to „On“ and the device begins to glow. After the application time the switch has to be switched off manually. During the application you can read, watch TV or carry out any other activities. A constant staring into the light source must be avoided.

Since this device emits a harmless light, it is not necessary to set up for this application a protection goggle. If you not like to look in a colored light, it is possible to wear a protection goggle during the application.

If you want to cover the phototherapy instrument after use allow the instrument to cool down at least 15 minutes.

Important: In order to be protected against UV-rays, the light box must only be operated with its original protective filter in place.

Fastening of light therapy units on the wall

The light therapy units can be mounted on the wall.

On the back side of the top of these light therapy units there is an opening in the middle below the air vents

(see illustration).

At this opening the units can be fixed on the wall with a screw. Depending on the type of wall, an appropriate screw must be used, may be with pegs. E.g. a round head wood screw which is screwed into an expansion anchor. The screw and the anchor are not included in the scope of delivery. The screw head must not exceed a max. height of 6 mm and a max. diameter of 8 mm. The distance from the wall to the screw head can be max. 8 mm. The screw shaft may have a maximum diameter of 4 mm. The load on the screw and the condition of the wall, to which the unit has to be mounted

must have a capacity of at least 3 kg.

The **Collagen-Light therapy unit BeautyLite BL 70** has on the back a fixture to swing out. This has to be removed before the fastening. Therefore 4 screws must be removed completely. Afterwards the screws have to be screwed in for safekeeping.

The base of the **Collagen-Light therapy unit BeautyLite BL 110** must be unscrewed before the fastening. Therefore the right and the left thumb screw must be removed completely from the foot. So the foot can be removed completely from the upper part of the unit. Afterwards the screws have to be screwed in for safekeeping.

Frequency and time of application, side effects

The application in the practice should be carried out 3 – 5 times per week for 20 minutes by the BL 110 and 30 minutes by the BL 70 until the skin appearance has improved. Subsequently, as required, as soon as the symptoms appear again. At home, the treatment can take place daily.

By using phototherapy the following side effects have been observed: Headaches, nausea, irritability and burning of the eyes can occur rarely for a short period. An absolute contraindication that means circumstances that do not allow the use of light therapy units is not known so far.

Please take care in case of all types of eye diseases and the eye specialist should be consulted before. In case of patients being under treatment with lithium or medicine for example Fluoxetine, Imipramine et al., which increase the sensibility of the eyes on light a supervision of the light therapy by an eye specialist could be considered.

These possible side effects / contraindications listed above can be avoided, if a protection goggle is worn during the application. If there is as side effect an itchy skin rash further treatments should be made only in consultation with the dermatologist.

Cleaning, maintenance and recycling

The instruments do not need any special maintenance. The manufacturer does not require such preventive inspections by other persons. The device requires no calibration. The device is not repairable and contains no user serviceable parts. No modification of this equipment is allowed. If the ME equipment is in use, servicing and maintenance is not allowed. The mains have to be disconnected, before opening the devices. To clean the outer casing parts of the devices, we recommend a mild cleaning agent (as used for glass).

Packaging: DAVITA® takes part in the system „Grüner Punkt“ („Green Point“). All packaging parts are made out of recyclable materials and must be recycled in accordance of the guideless concerning each material (eco-friendly). The casing and the front of the cover are made of plastic material. Fluorescent lamps as well as instruments marked with the symbol of a crossed out dustbin may not be disposed via the regular garbage.

The electrical and electronic components must be recycled separately. Please get in touch with your local authority to find the right way to of disposal. If electrical appliance is disposed of in landfills or dumps, hazardous substances can leak into the groundwater and get into the food chain, damaging your health and well-being. For the return and disposal of batteries the battery law is valid. National wide there are re-sampling points. The batteries must be thrown into the green collection boxes.

Technical data

Instrument Type:	Collagen-Light therapy unit BeautyLite BL 70	Collagen-Light therapy unit BeautyLite BL 110
Article number:	20500	20600
Housing:	High quality plastic housing	
Housing measures:	H: 52, W: 30, D: 7 cm	H: 52, W: 30, D: 7 cm
Weight:	2,1 kg	2,85 kg
Lamps:	2 x 36 watt Special lamps with red fluorescent tubes	3 x 36 watt Special lamps with red fluorescent tubes
Light spectrum:	580 – 630 nm	
Average lifetime of the lamps via the specification of the producer:	approx. 1.000 hours	approx. 1.000 hours
Energy costs: (at an average price of 0,20 € / kWh):	< 0,02 € per hour	0,02 € per hour
Fuse:	1 x fuse, 250 V ~ 3.15A Φ 5 x 20 The fuse sits behind the front housing	
Electric ballast units:	Instant starting, flickerfree light, no stroboscopic effects. Buzz free operation. Automatic switch off at the end of lamp live time.	
Light emission:	The selection of lighting source and light filter reduces the energy of light emitted in the Infrared and UV range to a minimum level, so it is harmless to human skin and eyes.	
Storage:	5 degrees C up to 50 degrees C and a relativ humidity of max. 93%, air pressure 700 to 1060 hPa	
Operation in the altitude of max.:	2.000 meter	
Operation temperature, humidity and air pressure:	10° C to 25° C; 15% to 93% relative humidity, air pressure 700 to 1060 hPa	
Mains requirements:	230 - 240 V~ , 50 - 60 Hz	

Signs and symbols

The following symbols appear in this manual, on the device or on its accessories:

Caution: Consult accompanying documents

Refer to instruction manual / booklet

Keep dry

CE Mark

Class II equipment

Manufacturer

SN

Serial number

DISPOSAL: Do not dispose this product as unsorted municipal waste. Collection of such waste separately for special treatment is necessary.

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Technical changes can be applied at any time without prior notice. Rev. No.: 11071511

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DAVITA® Medizinische Produkte GmbH & Co. KG • Im Hammereisen 29 A • D-47559 Kranenburg