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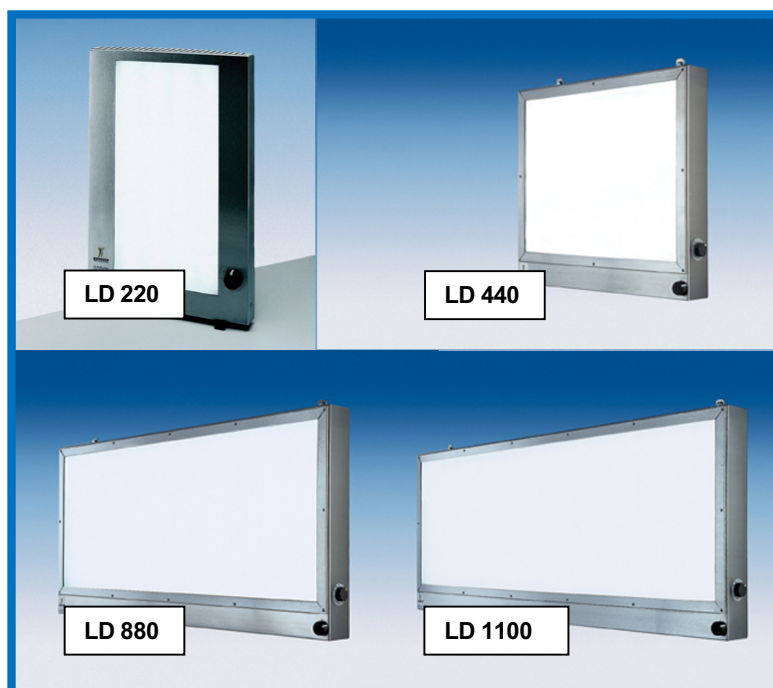
Dear customers,

We are pleased to welcome you as a new customer to our product. You have made the right decision. You have purchased a product which will improve your quality of life for many years.

Before you use the instrument, please read carefully this user manual.

## Medical explanation of the effect of light therapy

Today we know that sun light acts like a natural clock for many of our physiological processes. Bright light travels in the form of electric impulses via the eyes and the optical nerves to the nervous system's inner timer in our brain. This "time giver" processes the information about light and darkness and instructs other parts of the brain to rhythmically distribute chemical messengers. Bright light thus stimulates the production of chemical messengers which influences the mood and suppresses the production of the body's sleep inducing hormone.

## How to use the DAVITA® light therapy units

The light therapy boxes PhysioLight LD 220, LD 440, LD 880 and LD 1100 can be used against the following disease:

- ✚ Winter depression

## Some basic information before using the light therapy units

Before use, kindly check that the instrument is in functional order and that there is no sign of damage. The user must check that the equipment functions safely and see that it is in proper working condition before being used. Do not position the equipment so that it is difficult to operate the disconnection device. ☂ Phototherapy instruments must only be operated in dry and closed rooms that are free of any explosion risks. The device is not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide. The instruments must be standing free and should not be exposed to sun.

## How to carry out light therapy with the light therapy device PhysioLight LD 220

Place the instrument on an even surface and connect the plug with a wall socket. Sit down in front of the instrument and keep a distance of approx. 50 cm. The treatment time is at the above mentioned distance approx. 30 minutes. If you turn the control knob of the timer clockwise the light therapy device begins to glow. After the preset application time has expired, the LD 220 switches off automatically.

Keep your eyes open during the treatment. You are in the position to read, watch TV or carry out any other activities. After every 1 – 2 minutes you should look for a few seconds into the light source. A constant staring into the light source must be avoided. After the preset application time has run out, the instrument switches off automatically. It is also possible to switch off the instrument by turning the timer to zero

If you want to cover the phototherapy instrument after use allow the instrument to cool down at least 15 minutes.

Important: In order to be protected against UV-rays, the light box must only be operated with its original protective filter in place.

## How to carry out light therapy with the light therapy devices PhysioLight LD 440, LD 880 and LD 1100

One or two person(s) can be simultaneously treated with the LD 440; up to three persons can be simultaneously treated with the LD 880 and up to four persons with the LD 1100. There must be a distance of approx. 65 cm to the devices. The application time at this distance and with maximum luminous power is about 30 minutes.

By turning the control knob of the timer clockwise on the device's side, the light therapy unit begins to glow. First the selector switch on the instrument's front must be turned into one of the following positions:

Position 1: every 2nd lamp is on

Position 2: all lamps are on

Position 0: no lamp is on, even when the timer is on.

Keep your eyes open during the treatment. You are in the position to read, watch TV or carry out any other activities. After every 1 – 2 minutes you should look for a few seconds into the light source. A constant staring into the light source must be avoided. After the preset application time has run out, the instrument switches off automatically. It is also possible to switch off the instrument by turning the timer to zero

**Important:** In order to be protected against UV-rays, the light box must only be operated with its original protective filter in place.

## Operation of the light therapy units with the DAVITON®

The activated DAVITON® (not included) reminds you with an acoustical signal every 90 sec for 3 sec to look into the light source.

## Mounting the light therapy devices LD 440, LD 880 and LD 1100 on the wall

The light therapy devices Physiolight LD 440, the LD 880 and the LD 1100 can be fastened with screws and dowels on a flat wall (a distance of 505 mm for the LD 440, a distance of 1105 mm for the LD 880 and a distance of 1405 mm for the LD 1100).

A sufficient sustainability of the wall is necessary. You may have to use special fixings (for example if drywall walls are used). The installation has to be carried out by a qualified electrician. Under the device a free space of at least 20 cm must be left. Plug the female connector of the supplied power supply cable into the connector, located on the right bottom of the housing. Insert the male connector into a schuko plug socket.

## Frequency and time of treatment

If the treatment takes place in a clinic or in the doctor's surgery it should initially be carried out 5 times per week for a period of 2 – 3 weeks; thereafter, 3 times per week for a period of 2 – 3 weeks. After that, the treatment should take place as needed or when the symptoms occur again. At home, the treatment can take place daily, especially during autumn and winter time. The best treatment time for winter depression is the early morning or the time of the lowest mood.

## Side effects and contraindications

By using phototherapy the following side effects have been observed: Headaches, nausea, irritability and burning of the eyes can occur rarely for a short period. In case of porphyria, i.e. a metabolic disorder that is associated with a disturbance of the formation of red blood pigment heme, it is not allowed to treat with light therapy. Other contraindications, i.e., circumstances that prohibit the use of light therapy, are not known so far.

Please take care in case of all types of eye diseases and the eye specialist should be consulted before. In case of patients being under treatment with lithium or medicine for example Fluoxetine, Imipramine et al., which increase the sensibility of the eyes on light a supervision of the light therapy by an eye specialist could be considered.

## Cleaning, maintenance and recycling

The instruments do not need any special maintenance and requires no calibration. The manufacturer does not require such preventive inspections by other persons. Qualified companies / individuals can make repairs to the light therapy devices. The necessary technical documents and other information can be requested by DAVITA®. The mains have to be disconnected, before opening the devices. To clean the outer casing parts of the devices, we recommend a mild cleaning agent (as used for glass).

Packaging: DAVITA® takes part in the system „Grüner Punkt“ („Green Point“). All packaging parts are made out of recyclable materials and must be recycled in accordance of the guideless concerning each material (eco-friendly). The PhysioLight LD 220, LD 440, LD 880 and LD 1100 have a casing out of stainless steel.

Fluorescent lamps as well as instruments marked with the symbol of a crossed out dustbin may not be disposed via the regular garbage. The electrical and electronic components must be recycled separately. Please get in touch with your local authority to find the right way to of disposal. If electrical appliance is disposed of in landfills or dumps, hazardous substances

can leak into the groundwater and get into the food chain, damaging your health and well-being. For the return and disposal of batteries the battery law is valid. National wide there are re-sampling points. The batteries must be thrown into the green collection boxes.

## Technical data

Instrument Type:	Light therapy box PhysioLight LD 220	Light therapy box PhysioLight LD 440	Light therapy box PhysioLight LD 880	Ligth therapy box PhysioLight LD 1100
Article number:	10220	10440	10880	11100
EAN-Number:	4038962102209	4038962104401	4038962108805	4038962111003
Simultaneous treatment of:	1 person	Up to 2 persons	Up to 3 persons	Up to 4 persons
Medical product regarding 93/42/EWG	Yes	Yes	Yes	Yes
Housing:	High quality housing out of stainless steel			
Housing measures:	H: 60, W: 40, D: 8 cm	H: 66, W: 70, D: 10 cm	H: 66, W: 130, D: 10 cm	H: 66, W: 160, D: 10 cm
Weight:	7,5 kg	20 kg	36 kg	41 kg
Lamps:	4 x 55 watts	8 x 55 watts	16 x 55 watts	20 x 55 watts
Lighting power, distance and therapy length:	10.000 Lux at 50 cm distance for 30 minutes	10.000 Lux at 65 cm distance for 30 minutes	10.000 Lux at 65 cm distance for 30 minutes	10.000 Lux at 65 cm distance for 30 minutes
	5.000 Lux at 75 cm distance for 60 minutes	5.000 Lux at 100 cm distance for 60 minutes	5.000 Lux at 100 cm distance for 60 minutes	5.000 Lux at 100 cm distance for 60 minutes
	2.500 Lux at 110 cm distance for 2 hours	2.500 Lux at 150 cm for 2 hours	2.500 Lux at 150 cm for 2 hours	2.500 Lux at 150 cm for 2 hours
Average lifetime of the lamps via the specification of the producer:	approx. 10.000 hours	approx. 10.000 hours	approx. 10.000 hours	approx. 10.000 hours
Energy costs: (at an average price of 0,20 € / kWh):	0,04 € pro Stunde	< 0,07 € pro Stunde	< 0,15 € pro Stunde	< 0,19 € pro Stunde
Fuse:	2 x Fuses, 5 x 20 / T 2A H 250 V	2 x Fuses, 5 x 20 / T 3,15 A H 250 V	2 x Fuses, 5 x 20 / T 5 A H 250 V	2 x Fuses, 5 x 20 / T 6,3 A H 250 V
Electric ballast units:	Instant starting, flickerfree light, no stroboscopic effects. Buzz free operation. Automatic switch off at the end of lamp live time.			
Light emission:	The selection of lighting source and light filter reduces the energy of light emitted in the Infrared and UV range to a minimum level, so it is harmless to human skin and eyes.			
Storage:	5 degrees C up to 50 degrees C and a relativ humidity of max. 93%			
Mains requirements:	230 - 240 V~, 50 - 60 Hz			
Product Classification:	1. Class I equipment 2. No applied part 3. Not category AP / APG equipment 4. Mode of operation: Continuous operation			

## Signs and symbols

The following symbols appear in this manual, on the device or on its accessories:


Caution: Consult accompanying documents

Refer to instruction manual / booklet

Keep dry

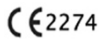CE Mark: conforms to essential requirements of the Medical Device Directive 93/42/EEC.

Class I equipment

Manufacturer

**SN**

Serial number

DISPOSAL: Do not dispose this product as unsorted municipal waste. Collection of such waste separately for special treatment is necessary.

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Technical changes can be applied at any time without prior notice. Rev. No.: 11101711

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DAVITA® Medizinische Produkte GmbH & Co. KG • Im Hammereisen 29A • D-47559 Kranenburg