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User Manual

Mode d'emploi

Magnetfeldtherapie-System Magneto Fit MF 400

Magnetic field therapy-System Magneto Fit MF 400

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Magnet field therapy-System Magneto Fit MF 400

Table of Contents-

Summary	6
Before operation	6
Safety instructions	6
Contraindications	7
Operation of Magneto Fit MF 400	7
Cleaning, maintenance and disposal	8
Guarantee	8
Technical specifications, scope of delivery and accessory	8

Before operation

Dear customer,

We enjoy that you bought a product from DAVITA® and congratulate you to your choice.

 Before you use the instrument, please read carefully this user manual.

In this user manual particularly the functions of the Magnetic field system MF 400 are described in detail.

Safety instructions

- ✚ Place the therapy unit on a compact and plane surface. In the mat and wristband are cables and coils among others, so an excessive, incorrect buckling has to be avoided.
- ✚ The mat is connected via a cable to the control unit. Pulling on the cable has to be avoided. The cable has a strain relief but it is possible to pull the cable from the mat.
- ✚ Check by mains operation if the power supply has the necessary 100 – 240 V by 50 – 60 Hz.
- ✚ The control unit should only be used with the wall power supply which belongs to it. The operation unit should only be connected to the appropriate treatment mat.
- ✚ This therapy unit should only be used inside. Avoid rain or much humidity.
- ✚ Do not place a case with liquid on or beside the unit.
- ✚ The magnetic field therapy has not the same effect for all people. Therefore every user should determinate the appropriate duration of the application by themselves. Many people use it approx. 30 minutes. Information for adjusting the necessary frequencies you can find for example in the provided book "Magnetfeldtherapie" from Dr. med. Gebhard Gehring, on page 35 and 36.
- ✚ Don't buckle the cables too much and do not expose the cables to heat or chemicals. If a cable is defect, pull it out of the unit and let it check by authorized qualified personnel.
- ✚ Do not open the unit. All service has to be done by authorized qualified personnel.
- ✚ Take off all metallic objects (jewels, watches etc.) before doing the treatment.


Contraindications

Warning: The magnetic field therapy is not allowed to use by patients:

- ✚ with a pacemaker (ad hoc decision) • with an ICD (implantated cardiofibrillator)
- ✚ with cancer (for cases in which the illness can be dated back for less than 5 years) • in the childbearing period
- ✚ with bacterial, inflammable athropathy•
- ✚ with massive, feverish infectious disease

In these cases please contact your physician .

Operation of Magneto Fit MF 400

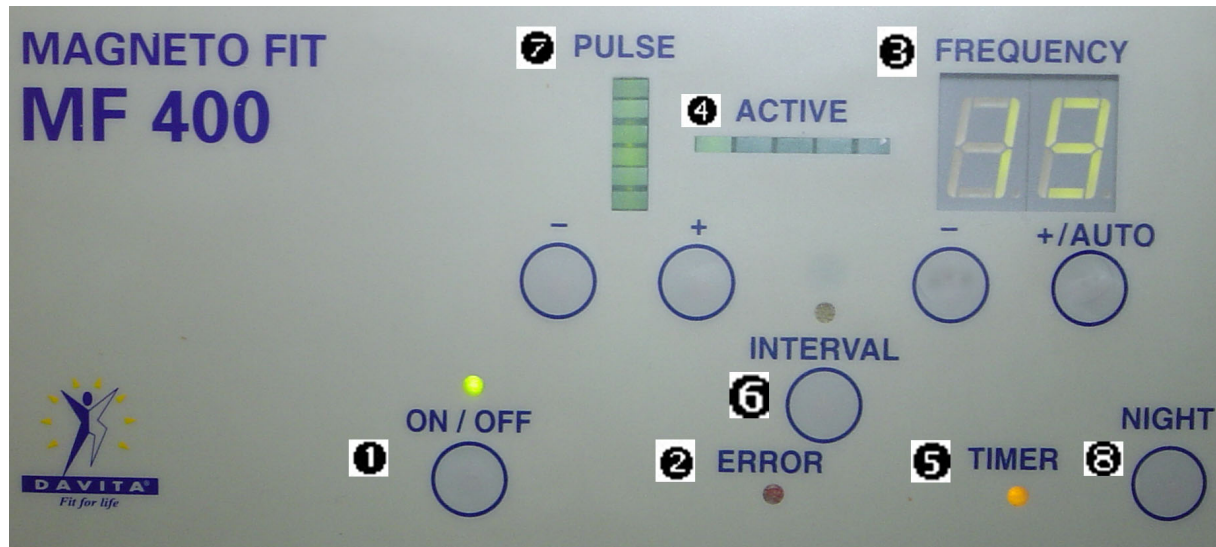Circuit points: ① Wall power supply, ② Magnetic field therapy mat

- ① **ON/OFF:** By pushing the button the unit can be switched on or off. After switching on the unit the ON/OFF-LED begins to shine.
- ② **ERROR:** This LED shines, if the mat is not or not correct connected with the control unit. The unit can not be used.
- ③ **FREQUENCY:** By pushing the „+“ button the frequency can be increased up to 19 Hertz. By pressing the „-“ button the frequency can be decreased until 0 Hertz. If 0 Hertz is reached the control unit is deactivated.
FREQUENZROTATION: If the „+“ button will be pushed a little longer at 19 hertz the frequency rotation function will be activated. Here the single frequencies are rotating active for 5 seconds. Many users like this function very much. By pushing the „-“ button this modus ca be left.
- ④ **ACTIVE:** If the control unit is operating, the LEDs of the display will be powered one after another (Running light). The speed of the running light is independent from the adjustment of the frequency.
- ⑤ **TIMER:** If the frequency is adjusted higher than 4 hertz the timer switches on and the LED begins to shine. The Timer will be on for 30 minutes. After this time the units switch off. The frequency will be adjusted to 0 Hertz.
- ⑥ **INTERVAL:** By pushing the button the LED will be activated. There is a change in the function of the unit. The control unit will be switched on for 6 seconds and after that switched off for 3 seconds (Rest) After that the process starts again. By pressing the button again the function can be deleted. Some scientists have the opinion, that the body can assimilate the energy of the magnetic fields so optimal.
- ⑦ **PULSE:** By pushing the „+“ or „-“ button the intensity will be increased or decreased. Here exists the possibility to variate the magnetic active time between 15% and 90%. The application time can be considerably shortened. Every

illuminated line in the LED means an increase of 15%.

- ③ **NIGHT:** By pushing the button all displays will be cleared. Only the small LED in the frequency display is shining. The function of the unit will be unaffected. By pressing any button the displays will be switched on again.

Storage of the operation parameters:

By pushing the ON/OFF-button the intensity, frequency and the frequency rotation modus will be stored. If the unit will be switched on again these parameters will be imported and the unit works with these parameters. So it is possible to switch on the unit with an usual clock timer and the unit starts with the stored parameter.

Cleaning, maintenance and disposal

A special maintenance is not necessary for this device. Before cleaning the unit has to be switched off and all plugs have to be pulled out. For cleaning the unit a soft, wet rag and a mild household cleaning agent can be used. The cleaning agent should not be biting or rubbing. The magnetic therapy mat is carefree.

Packaging: DAVITA® joint to the system „Grüner Punkt“. All parts of the packaging are recycable materials and environmentally acceptable to dispose. Batteries and electronic units which are signed with the crossed wheelie bin symbol are not allowed to throw into the domestic waste. There is the commitment to dispose such electro- and electronic „Altgeräte“ separately. The communal public authority which is responsible for the disposal gives information regarding the possibilities of a regulated disposal.

Note STK: In accordance to the “Medizinprodukte-Betreiberverordnung (MP BetreibV) § 6 (Sicherheitstechnische Kontrollen)” the operator is obliged to let carry out a safety technical control every 2 years. We advise to carry out the controls by DAVITA® or from a special authorized qualified personnel. The unit which should be maintained can be dispatched directly to the producer or to the distributor, where it has been sold.

Technical specifications, scope of delivery and accessory

Wall power supply

Model: UE 15WCP – 120125SPA
Input: 100 - 240 V AC 50-60 Hz
Output: 12 V DC 1,25 A max.

Control unit

Type: MF 400, CE 2004
Power consumption: 12 V DC / 8,5 VA
Output: 12 V / 6,5 W
Frequency: 1 – 19 Hz
Magnetic induction: 0,21 mT max.
Automatic switch off:
From 1 – 4 Hz no switch off
From 5 – 19 Hz switch off after approx. 30 minutes

Applicator

Mat: 190 x 85 cm
Power consumption:
12 V pulsed 0,7 A
Field intensity: 0,21 mT max.

Weight of the complete system: approx. 4,2 kg

Scope of delivery: Control unit, mat, wall power supply and detail user manual.

Magneto Fit MF 400: Art. No.: 60010, EAN-No.: 4038962600101

Accessory:

Wristband out of cotton, 80 cm long and 20 cm wide, Art.No.: 62001

Treatment mat 60 cm long and 50 cm wide, Art.No.: 62006

The magnetic field therapy system Magneto Fit MF 400 is a proofed medical product in accordance to DIN EN 60601-1-2 – Made in Germany. That guarantees a high quality and the fulfillment of all the relevant safety standards. Manufacturer: Hane-Produktmanagement GmbH, D-86156 Augsburg

On general the magnetic field therapy is approved and used from the scientific medicine. We are obliged by law to point out, that there are existing alternative opinions regarding the effect, therapy success and therapy improvement in the scientific medicine.

C E 2004

Subject to technical changes. Rev. No. 11071511

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