Lichttherapiegeräte

Light therapy boxes

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Lichttherapie apparaten

Valoterapialaitteet

Ljusterapilampor

VILUX® 70 - VILUX® 100, VITA^{bright}® 70 - VITA^{bright}® 100

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User manual

Mode d'emploi

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Bruksanvisning

Deutsch.....	2
English.....	6
Français.....	10
Nederlands.....	14
Suomi.....	18
Svenska.....	22

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Dear customers,

We are pleased to welcome you as a new customer to our product. You have made the right decision. You have purchased a product which will improve your quality of life for many years.

Medical explanation of the effect of light therapy

Today we know that sun light acts like a natural clock for many of our physiological processes. Bright light travels in the form of electric impulses via the eyes and the optical nerves to the nervous system's inner timer in our brain. This "time giver" processes the information about light and darkness and instructs other parts of the brain to rhythmically distribute chemical messengers. Bright light thus stimulates the production of chemical messengers which influences the mood and suppresses the production of the body's sleep inducing hormone.

How to use the DAVITA® light therapy units

The light therapy boxes can be used against the following disease:

- ✚ winter depression

Some basic information before using the light therapy units

Before use, kindly check that the instrument is in functional order and that there is no sign of damage. Phototherapy instruments must only be operated in dry and closed rooms that are free of any explosion risks. The device is not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide. The instruments must be standing free and should not be exposed to sun.

How to carry out light therapy

Place the instrument on an even surface and connect the plug with a wall socket. The light therapy boxes have on the back a fixture to swing-out, so that the units have a save stand. Sit down in front of the instrument and keep a distance of approx. 55 cm in case of the VILUX® 70 / VITAbright® 70 and in a distance of approx. 45 cm in case of the VILUX® 100 / VITAbright® 100. The treatment time with the VILUX® 70 / VITAbright® 70 is at the above mentioned distance approx. 2 hours. The treatment time with the VILUX® 100 / VITAbright® 100 is at the above mentioned distance approx. 1 hour.

Move the tilting switch to „On“ and the device begins to glow. After the application time switch off the device.

Keep your eyes open during the treatment. You are in the position to read, watch TV or carry out any other activities. After every 1 – 2 minutes you should look for a few seconds into the light source. A constant staring into the light source must be avoided. If you want to cover the phototherapy instrument after use allow the instrument to cool down at least 15 minutes.

Important: In order to be protected against UV-rays, the light box must only be operated with its original protective filter in place.

Fastening of light therapy units on the wall

The light therapy units can be mounted on the wall.

On the back side of the top of these light therapy units there is an opening in the middle below the air vents (see illustration).

At this opening the units can be fixed on the wall with a screw. Depending on the type of wall, an appropriate screw must be used, may be with pegs. E.g. a round head wood screw which is screwed into an expansion anchor. The screw and the anchor are not included in the scope of delivery.

The screw head must not exceed a max. height of 6 mm and a max. diameter of 8 mm. The distance from the wall to the screw head can be max. 8 mm. The screw shaft may have a maximum diameter of 4 mm. The load on the screw and the condition of the wall, to which the unit has to be

mounted must have a capacity of at least 3 kg.

The light therapy units have on the back a fixture to swing out. This has to be removed before the fastening. Therefore 4 screws must be removed completely. Afterwards the screws have to be screwed in for safekeeping.

Operation of the light therapy units with the DAVITON®

The activated DAVITON® (not included in the scope of delivery) reminds you with an acoustic signal every 90 seconds to look for 3 seconds direct into the light source.

Frequency and time of treatment

At home, the treatment can take place daily, especially during autumn and winter time. The best treatment time for winter depression is the early morning or the time of the lowest mood.

Side effects and contraindications

By using phototherapy the following side effects have been observed: Headaches, nausea, irritability and burning of the eyes can occur rarely for a short period. In case of porphyria, i.e. a metabolic disorder that is associated with a disturbance of the formation of red blood pigment heme, it is not allowed to treat with light therapy. Other contraindications that means circumstances that do not allow the use of light therapy units are not known so far.

Please take care in case of all types of eye diseases and the eye specialist should be consulted before. In case of patients being under treatment with lithium or medicine for example Fluoxetine, Imipramine et al., which increase the sensibility of the eyes on light a supervision of the light therapy by an eye specialist could be considered.

Cleaning, maintenance and recycling

The instruments do not need any special maintenance and requires no calibration. The manufacturer does not require such preventive inspections by other persons. Qualified companies / individuals can make repairs to the light therapy devices. The necessary technical documents and other information can be requested by DAVITA®.

The mains have to be disconnected, before opening the devices. To clean the outer casing parts of the devices, we recommend a mild cleaning agent (as used for glass).

Packaging: DAVITA® takes part in the system „Grüner Punkt“ („Green Point“). All packaging parts are made out of recyclable materials and must be recycled in accordance of the guideless concerning each material (eco-friendly).

The casing and the front of the cover are made of plastic material. Fluorescent lamps as well as instruments marked with the symbol of a crossed out dustbin may not be disposed via the regular garbage.

The electrical and electronic components must be recycled separately. Please get in touch with your local authority to find the right way to of disposal. If electrical appliance is disposed of in landfills or dumps, hazardous substances can leak into the groundwater and get into the food chain, damaging your health and well-being. For the return and disposal of batteries the battery law is valid. National wide there are re-sampling points. The batteries must be thrown into the green collection boxes.

Technical data


Instrument Type:	Light therapy box VILUX® 70 / VITAbright® 70	Light therapy box VILUX® 100 / VITAbright® 100
Article number:	13070	13100
Medical product regarding 93/42/EWG:	Yes	Yes
Housing:	High quality plastic housing	
Housing measures:	H: 52, W: 30, D: 7 cm	H: 52, W: 30, D: 7 cm
Weight:	2,2 kg	2,2 kg
Lamps:	2 x 36 watt	3 x 36 watt
Lighting power, distance and therapy length: Important: Please notice that the light therapy can only successfully used in a distance of at least 45 - 55 cm!	10.000 Lux at 20 cm distance for 30 minutes	10.000 Lux at 30 cm distance for 30 minutes
	5.000 Lux at 35 cm distance for 60 minutes	5.000 Lux at 45 cm distance for 60 minutes
	2.500 Lux at 55 cm distance for 2 hours	2.500 Lux at 70 cm distance for 2 hours
Average lifetime of the lamps via the specification of the producer:	approx. 10.000 hours	approx. 10.000 hours
Energy costs: (at an average price of 0,17 € / kWh):	< 0,02 € per hour	< 0,02 € per hour
Fuse:	1 x fuse, 250 V ~ 3.15A Φ 5 x 20 The fuse sits behind the front housing	
Electric ballast units:	Instant starting, flickerfree light, no stroboscopic effects. Buzz free operation. Automatic switch off at the end of lamp live time.	
Light emission:	The selection of lighting source and light filter reduces the energy of light emitted in the Infrared and UV range to a minimum level, so it is harmless to human skin and eyes.	
Storage:	5 degrees C up to 50 degrees C and a relativ humidity of max. 90%	
Operation in the altitude of max.:	2.000 meter	
Operation temperature and humidity::	10° C to 25° C; 10% to 95% relative humidity	
Mains requirements:	230 - 240 V~ , 50 - 60 Hz	
Product Classification	<ol style="list-style-type: none"> 1. Class II equipment 2. No applied part 3. Protection against ingress of water: IP 20 4. Not category AP / APG equipment 5. Mode of operation: Continuous operation 	

Signs and symbols

The following symbols appear in this manual, on the device or on it's accessories:

Caution: Consult accompanying documents

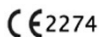CE Mark: conforms to essential requirements of the Medical Device Directive 93/42/EEC.

Class II equipment

Manufacturer

SN

Serial number

DISPOSAL: Do not dispose this product as unsorted municipal waste. Collection of such waste separately for special treatment is necessary.

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Technical changes can be applied at any time without prior notice. Rev. No.: 11081511

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DAVITA® Medizinische Produkte GmbH & Co. KG • Im Hammereisen 29 A • D-47559 Kranenburg